

# BIOPSY



Two trigger buttons: one on the side and one on the rear



Somacore Verify is tailored to meet a wide range of clinical requirements during biopsy procedures. The advanced delay mode enables physicians to reposition the biopsy cannula before cutting the tissue.

Somacore Verify		
Order No.	Size	Box Qty.
900380-US	14 G x 10 cm	5 units
900390-US	16 G x 10 cm	5 units


Puncture Sheath	
Order No.	Box Qty.
170531-US	10 units
170541-US	10 units

# Somacore Verify

For Soft Tissue Biopsy



Atraumatic needle tip and echogenic marking for precise Ultrasound positioning

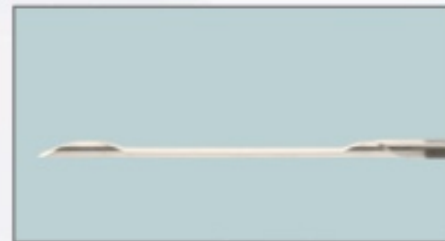

Large 18 mm sample notch

## FEATURES:

- Three unique firing positions
  - No throw mode
  - Delay mode
  - Automatic firing mode
- Double Diamond Arrow tip point designed to penetrate hard lesions
- Echogenically 'etched' to ensure optimal visibility under Ultrasound
- Reliable quality of biopsy samples by a large sample notch and sharp cutting cannula
- Innovative delay mode to allow tissue to settle in the sample notch and reposition the cannula before cutting if necessary
- Two trigger buttons allow for different biopsy handling techniques according to each clinical case